

NatriFlo® HD-Q Membrane Adsorber:

A new single-use high-capacity polishing tool with excellent salt tolerance and process robustness for mAb purification

NATRIX SEPARATIONS

Burlington, ON, Canada

Dharmesh Kanani, Ph.D.

Renaud Jacquemart, Ph.D.

Navneet Sidhu, Annabel Shang,

Mochao Zhao, James Stout, Ph.D.

For more information:

Dharmesh Kanani

Natrix Separations

kanani@natrixseparations.com

ABSTRACT


Strong anion exchange (Q) chromatography has become an industry standard in mAb production. It is a proven technology to remove DNA, viruses, endotoxins and acidic host cell proteins from process feed streams in flowthrough mode.

Recent trends show an increasing interest in downstream single-use technologies and flexible biomanufacturing due to advancement in cell culture technology and emergence of biosimilars. Traditional chromatography columns are slow, often oversized and not suitable for flexible biomanufacturing. Conventional membrane adsorbers cannot provide sufficient process robustness due to low binding capacity of membranes. These factors impose challenges on the design of purification schemes for manufacturing of biotherapeutics. However, the NatriFlo HD-Q membrane adsorbers overcome these limitations by combining high binding capacity and high flow rates to deliver best-in-class performance in single-use plug & flow format.

OBJECTIVES: This study compares dynamic BSA binding capacity of NatriFlo HD-Q with competing technologies to demonstrate performance in challenging process conditions (salt & phosphate background). It also demonstrates the effectiveness of NatriFlo HD-Q in removing process impurities (HCP, DNA & viruses) from industrial feed streams.

INTRODUCTION

Three Chromatography Steps in mAb Production Platform

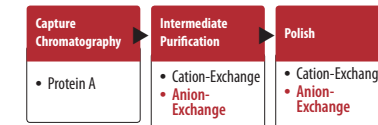

Anion Exchange (AEX): Process Requirements

- Better & robust separation performance with increased throughput and capacity in flowthrough mode
- Flexibility in process design (such as positioning of AEX step in the platform, salt tolerance and buffer compatibility)
- Risk-free, scalable & easy operation in single-use format

AEX: Pros/Cons of Existing Technologies

Column chromatography:

- Robust separation performance due to high media binding capacity but often ends up being oversized due to throughput limitation
- Relatively complex operation & not amicable to single-use format

Membrane chromatography:

- High throughput & good performance only when feed is relatively clean due to limited media binding capacity
- Lower process robustness and less design flexibility

NatriFlo HD-Q Membrane Adsorber Overview

- NatriFlo HD-Q is a new high capacity, high throughput strong anion exchange membrane adsorber designed to purify biomolecules such as mAbs in flow-through mode
- Features advanced membrane material with a 3D macroporous hydrogel structure that provides High Density of binding sites and rapid mass transfer
- NatriFlo HD-Q Family offers scalable polishing from R&D to cGMP manufacturing

CONCLUSIONS

NatriFlo HD-Q Benefits: Reduced Risk, Greater Flexibility & Improved Economy

- Improved salt tolerance and wider buffer compatibility**
 - Superior dynamic BSA binding capacity in the presence of salt or phosphate background
- High loading capacity without compromising process robustness**
 - Excellent virus and DNA clearance from partially purified mAb sample for load as high as 10 kg/L over a broad design space
 - Excellent HCP clearance with no sign of increasing breakthrough even at 20 kg/L load
- Simple scale-up**
 - Excellent & consistent viral and HCP clearance from lab to process scale
- Truly single-use plug & play format**
 - No column packing, qualification, cleaning or validation required, thus significantly reducing capital expenditures and operating expenses
 - No risk of cross-contamination


COMPETITIVE EVALUATION

Experimental


DEVICE	MATRIX	CHEMISTRY	MEDIA VOLUME (mL)	FLOW RATE (mL/min)
NatriFlo HD-Q Recon Mini	Porous Polyacrylamide	Q	0.20	2.0
Membrane 1-Q	Modified hydrophilic polyethersulfone	Q	0.18	1.8
Membrane 2-Q	Stabilized reinforced cellulose		0.08	0.8
Resin 1-Q	Highly cross-linked agarose		1.00	1.0
Resin 2-Q	6% Highly cross-linked agarose		1.00	1.0
Membrane 3-PA	Stabilized reinforced cellulose	Primary Amine	0.08	0.8
Membrane 4-PA	Ultra high MW polyethylene		0.08	0.8
Resin 3-MM	Highly cross-linked agarose		Mixed Mode	1.00

- All the devices were run as per the instructions provided by the vendors at recommended flow rate (i.e. membrane adsorbers were run at 10 MV/min whereas 1 mL prepacked columns were run at 1 CV/min)
- Feed sample:
 - 1 g/L BSA in 25 mM Tris + NaCl buffer, pH 8.0 (salt tolerance study)
 - 1 g/L BSA in phosphate buffer, pH 8.0 (phosphate tolerance study)
 - 0.1 g/L Herring testes DNA in 25 mM tris buffer, pH 8.0 (DNA binding study)


10% BT BSA Binding Capacity: Salt Tolerance

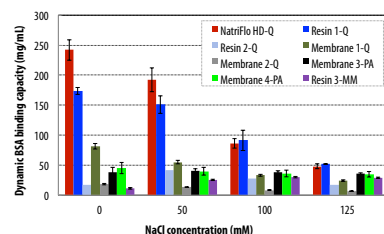

10% BT BSA Binding Capacity: Phosphate Tolerance

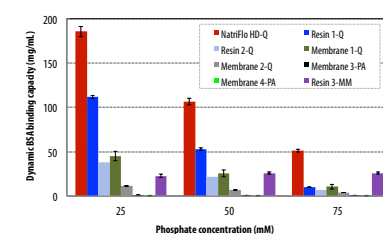

10% BT DNA Binding Capacity

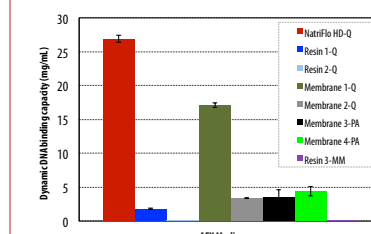


- The superior dynamic binding capacity of NatriFlo HD-Q membrane adsorber over a wide range of operating conditions combined with impressive throughput enable simple, robust, low cost polish operations in a single-use format.
- Resin 1-Q achieved adequate dynamic BSA binding capacity but very low dynamic DNA binding capacity at 1 minute of residence time.
- In most cases, column chromatography ends up being oversized due to throughput limitation. For example, a typical column would need > 66 hrs just to load the feed sample having 5 g/L titer to 10 kg/L capacity at 2 minutes of residence time. On the other hand, NatriFlo HD-Q would need only 3.3 hours to load the same feed at 10 MV/min flow rate (or 1.3 hours at 25 MV/min flow rate).
- Conventional membrane adsorbers lack the process robustness & design flexibility due to poor binding capacity of media.

COLLABORATOR'S DATA

Virus/DNA Clearance & Process Robustness

Robust, Load Independent MVM Clearance

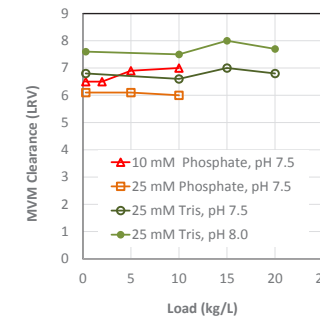

Load independent MVM clearance over wide design space


Load: 0.25 – 20 kg/L membrane (0.25 kg/L typical resin load)

Phosphate vs. Tris


Sample: 10 g/L mAb, 10 mS/cm

Residence time: 6 seconds

Similar and better MVM Clearance at pH 6.5-7.0

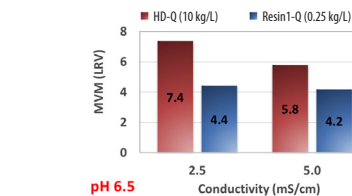


HD-Q gives higher MVM Clearance than state of the art Q resin
40x smaller disposable device has similar processing times, gives equal or superior results

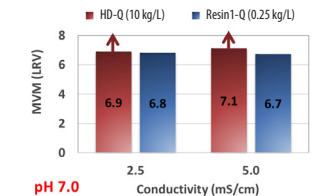


	HD-Q	Resin1-Q
Load	10 kg/L	0.25 kg/L
Residence time	0.1 minute	3 minutes
Sample: 10 g/L mAb in 10 mM Phosphate + NaCl		
Note: Up arrow (↑) indicates limit of detection		

HCP Reduction Performance at Lab Scale for 10 kg/L Load

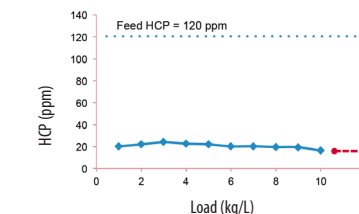

Sample: 8 g/L Protein A purified mAb in 20 mM phosphate + 100 mM NaCl, pH 7
Flow rate: 10 MV/min (residence time – 6 seconds)

Linear Scalability

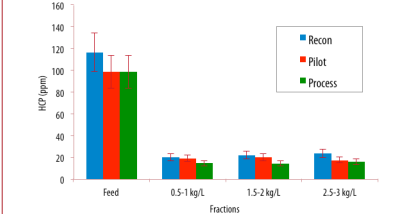

Trial at collaborator (Courtesy of Patheon)

- Excellent HCP clearance with no sign of increasing breakthrough even at 10 kg/L load (lab scale)
- Excellent & consistent HCP clearance from lab to process scale, even in phosphate buffer (Pilot & Process devices were tested only up to 3 kg/L load due to limited supply of protein A purified mAb)
- Consistent DNA removal (>2.8 LRV) across all scales for 3 kg/L load (data not shown)